

Moving & Handling

Transfer Boards	752 - 755
Handling Belts	756
Glide Sheets	757 - 758
Patient Turners	759 - 761
Turntables	762 - 763
SuperPole™	764
Bed Assists	765 - 766
Seat Risers	767
Patient Transfer	768



BeasyTrans™ Easy Transfer Systems

This product range utilises the innovative idea of a rotating seat that tracks along a rigid plastic board. This reduces the risk of shear stress and friction associated with patient transfer because the seat moves with the user. The user glides laterally on a round seat that rotates through 360°, meaning that it is easy to use and highly practical for the user. The BeasyTrans Transfer Systems are made from exceptionally strong plastic, which means that they can support a user up to 180kg (28¼st). They are also light enough to be easily transportable. The Beasy transfer system is available in three styles.

Maximum user weight



1 BeasyTrans™

This transfer system is ideal for car to wheelchair transfers. Length 1020mm (40"). Seat diameter 300mm (12"). Weight 2.9kg. **081288323**

2 BeasyGlyder™

Especially designed for transfer past wheelchair arms and similar obstacles. Length 810mm (32"). Seat diameter 250mm (10"). Weight 2.2kg. **081288331**

3 BeasyII™

Ideal for short transfer distances, such as between bed and wheelchair. Length 700mm (27½"). Seat diameter 300mm (12"). Weight 2kg. **081288349**

1 Days Curved Transfer Board

The Days Curved Transfer Board is designed to enable people with limited strength and mobility to assist with functional transfers. The board is ideal for car, bed, chair and toilet transfers, facilitating independence with safety. The arc shape not only enables a variety of positioning options, but also allows for the natural sweeping pattern of an assisted transfer. For increased stability and security the Days Curved Transfer Board has two non-slip pads on the base of the board. Wipes clean easily for infection control purposes. Length 710mm (27½"). Width 370mm (14½"). Weight 1.95kg.

081306554

Maximum user weight



2 Q2 Banana Board

The Q2 Banana Board is a truly universal transfer board. Its unique curved shape and maximum working load of 285kg enable it to be used in virtually all applications. Ideal for use in both assisted and non-assisted transfers, it can also be used for face-to-face transfers within restricted areas. The inset anti-slip pads on the underside prevent pad detachment and afford maximum grip on the two transfer surfaces. It is lightweight and is easily wiped clean, using either household disinfectants or alcohol-based solutions. Length 750mm (29½"). Width 338mm (15¼"). Weight 1.6kg

091120534

Maximum user weight




B

3 Duo Transfer Slide Deluxe

The unique design of the Duo Transfer Slide Deluxe allows for an angled transfer as well as a straight transfer from one seated position to another. The handgrips at either end give the user a means of pulling themselves across the smooth surface whilst ensuring the slide is securely located during the transfer. Length 750mm (29½"). Width at centre 220mm (8½").

091120468

Maximum user weight

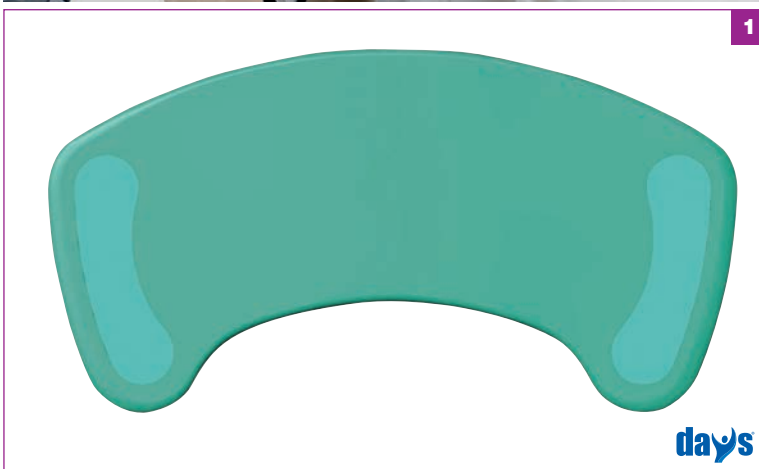

B BARIATRIC
These products are designed for bariatric users

1



days

1 Days Wooden Transfer Boards

These durable transfer boards are ideal for bed, wheelchair and car transfers. Varnished to facilitate a friction free transfer, reducing strain and eliminating the need to lift patients. The chamfered ends ensure secure positioning during transfers. The longer board has hand holes cut into both ends for easier grip.

	Length	Weight
091120633	600mm (24")	1.1kg
091120724	760mm (30")	1.38kg

Maximum user weight



2 Banana Board Lateral Folding Board

The Banana Board Lateral Folding Board has been specifically designed to enhance the versatility of Lateral Transfers. Its unique design enables the board to either match the contours of a profiling bed or be used completely flat as in a normal lateral transfer. The Banana Board's ability to be able to transfer directly across to a profiled bed is particularly useful where a patient is either afraid or unable to lie completely flat during a lateral transfer. Lightweight with six ergonomically positioned handles, the Banana Board is easily placed and removed in most transfer situations. The handholds, when the board is folded in half, can be used to hang the board up, thereby ensuring that the board takes up the minimum of space during storage. Recommended maximum gap width 150mm (6") Length 1500mm (59") Width 740mm (29").

091188085

1



days

1



days

2



1 Butterfly Board and Wing

A long, wide and stable transfer board with an optional GlideWing for additional slide function. The Butterfly board low friction upper surface allows the user to transfer independently or with assistance. Suitable for transferring between bed and wheelchair / shower chair, between wheelchair / chair / toilet chair or to and from the car. Combining the Butterfly board with the GlideWing provides even lower friction and is especially recommended when the transfer involves bare skin, for instance to and from a shower chair or toilet. The Butterfly board has cut outs for positioning against the wheel of a wheelchair and non-slip pads on the underside allowing for good positioning and preventing the board from sliding out of place. Washable at 75°C. Length 660mm (26"). Width 320mm (12½"). Weight 1.35kg.

091120914 Butterfly Board

091120922 GlideWing

Maximum user weight



2 Transfer Board Glide Sheet

Made from low friction material, this glide sheet is especially designed for use with transfer boards. It features a non-slip flap to prevent the glide sheet moving away from the transfer board whilst the transfer is in progress. Size 300mm x 700mm (12 x 27 ½").

091121763

3 Dolphin Board

The Dolphin Board is designed for assisted and independent transfer from a wheelchair to the toilet. It has the additional benefit of remaining in position during toilet use. The wheelchair cut-outs allow transfer from either side of the toilet. Strong and lightweight, the board is rigid with no flex during transfer. Length 690mm (27"). Width at widest point 380mm (15"). Weight 1.65kg.

091120971

Maximum user weight



4 Commode Transfer Board

Uniquely shaped with strategically placed cut-out, this board does not need to be removed during an individual's toileting routine after the client transfers to the commode or toilet. The white plastic, high density board has two different sides; one is striate-sanded for wet, bare skin transfers, the other is buff-sanded for clothed transfers. The birch transfer board is available in two weight limits. Length 480 (19"). Width 430mm (17"). Thickness 13mm (½"). Weight 1.8kg.

081229954

Maximum user weight



1

MAXI



MAXI PLUS



2

MAXI



3



1 Handling Belt

These generous padded belts are comfortable for the patient and allow the user to obtain a good secure grip. Benefitting from a series of vertical and lateral handholds which give many options when transferring or walking with the patient. Available in two sizes to suit a range of patients from children through to larger adults. Machine washable up to 71°C.

Waist Size

091121458 Maxi 61-114cm (24-45")
091121466 Maxi Plus 91-152cm (36-60")

2 Deluxe Handling Belts

Similar to the Handling Belts above but with a non-slip inner lining to help prevent the belt sliding upwards or sideways during transfer. Machine washable up to 71°C.

Waist Size

091121474 Maxi 61-114cm (24-45")

3 Universal Handling Belt

The Universal Handling Belt is designed to accommodate waist sizes from 51cm (20") to 145cm (57") and features front and side loops that are maintained at any extended length. It is manufactured from material that helps to reduce cross infection and is impregnated with an anti-bacterial agent. It can be washed up to 75°C.

091122084



Days Easi-Mover Glide Sheets

Designed from lightweight fabric, the Easi-Mover Glide Sheets have ultra low friction surfaces allowing easy movement of the patient. The Easi-Mover range includes tubes and flat sheets. We recommend that Easi-Movers be used in accordance with existing moving and handling policies or with the advice of your moving and handling advisor.

1 Days Easi-Mover III Glide Sheets

Here the fabric is doubled and sewn together in a tube, ensuring a smooth glide between two surfaces, thus reducing skin friction for the patient and strain on the carer. The sheets are yellow with colour-coded reinforcement webbing to differentiate size. Washable to 65°C.

Size	Webbing Colour
091120252 580×570mm (23×22½")	Grey
091120260 900×900mm (35½×35½")	Red
091120278 1200×710mm (47×28")	Royal Blue
091120286 1200×930mm (47×36½")	Green
091120294 1450×710mm (57×28")	Black
091120302 1900×740mm (74×29")	White


2 Days Easi-Mover Flat Sheets

Made from the same yellow fabric as the Easi-Mover III (above) but flat rather than looped and with handles to aid grip. Again, the colour-coded reinforcement edging is used to differentiate size. Washable to 65°C.

Size	Webbing Colour
091120310 900×1400mm (35½×55")	White
091120328 750×2000mm (29½×78¾")	Green
091120344 1000×2000mm (39½×78¾")	Royal Blue
091120351 1500×2000mm (59×78¾")	Red



£ PRICES
For up-to-date prices please see our Website, alternatively contact your Sales Representative or Customer Services.

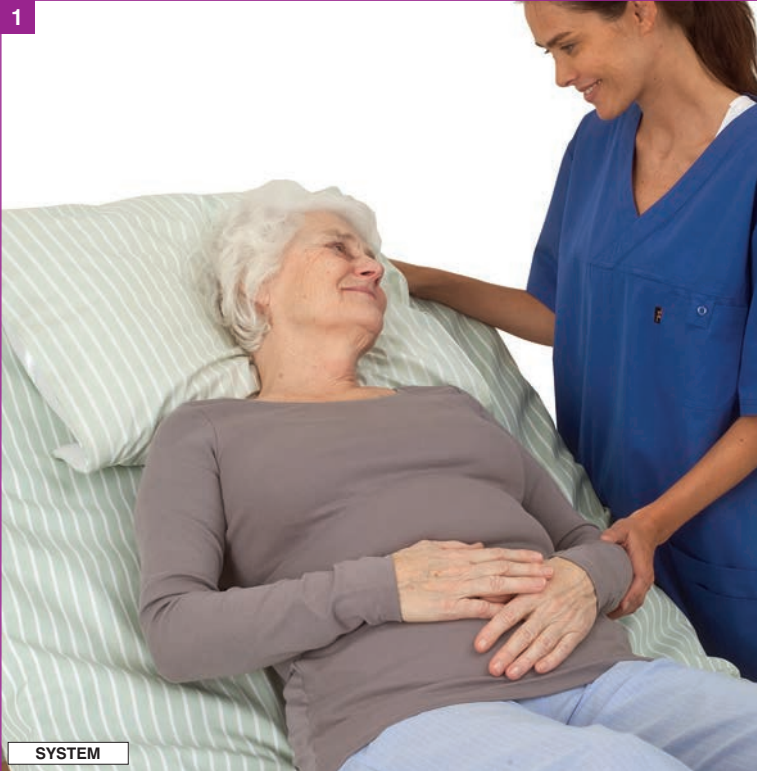

1 Etac® SatinSheet System

Etac® SatinSheet System is developed to cover many different levels of need. The system can be used for everyone from a user with slightly reduced/impaired mobility that may be able to turn independently in bed, to a user with high need of caregiver assistance.

SatinSheet 2Direction

SatinSheet 2D has a low-friction satin area in the centre of the sheet that allows movement from side-to-side thus facilitating turning in bed. The edges of the sheet are of cotton/polyester and will brake movement, allowing the user to sit on the edge of the bed without slipping.

SatinSheet, 2D Maxi (sliding surface 70×200cm)

SatinSheet 4Direction – combination for complete assisted movement SatinSheet 4D is a draw sheet with satin on the underside and specifically adapted to be used in combination with a SatinSheet 2D underneath. The combination of a SatinSheet 2D base sheet with a SatinSheet 4D creates a system that allows movement in all directions. This solution is suitable for a user with limited mobility.

SatinSheet 4D Midi (satin 148×140cm)

091530765 SatinSheet System - 2D Maxi + 4D Midi

Maximum user weight



2 Glide and Lock Sheets

Glide and lock sheets provide a low friction movement in one direction - into the chair or wheelchair or up the bed - but resist the natural tendency of the patient to slide down again. In normal use, the product allows repeated positioning whilst considerably reducing the frequency of the moves. This not only improves patient care but also decreases the handling risks faced by staff dealing with these types of manoeuvre on a regular basis. Glide and lock sheets are an invaluable resource in any department or homecare environment.

Features:

- Internal locking strips to improve client safety
- Anti-slip material on the base to prevent unwanted movement of the sheet
- Variety of sizes and finishes to match the product with the client's needs
- Washable at 72°C

	Type	Size
091120559	Wheelchair	430×450mm (17×18")
091120567	Chair	480×530mm (19×21")
091120583	Bed	650×800mm (26×32")



WHEELCHAIR



CHAIR



BED

1 Rota Stand Compact

The folding and lightweight nature of the Rotastand Compact is ideal for community assessments by therapists who often have to lift equipment in and out of their vehicles. Equally, a more independent lifestyle is encouraged for clients who wish to travel to friends and relatives. The Rotastand Compact fits easily into the boot of a car. Features include a height adjustable frame and shin pad to suit all clients; central wheels to allow easy positioning and a stabilizing bar to allow a one person transfer. Weight 10kg.

091544527

Maximum user weight



2 Rota Stand

Developed to allow an assisted transfer from one seated position to another. The frame, shin pad, and turntable combine to give optimum safety to both the carer and the user during the transfer. The frame and shin pad are height adjustable to suit all and a foot pedal break to secure the turntable. The central wheels allow the Rota Stand to be rolled up to the user even after their feet have been positioned on it. Handle height 790 to 1300mm (31 to 50"). Weight 15kg.

091120765

Maximum user weight



3 Rota Stand Solo

The Rota Stand Solo features a special mechanism in the base that allows the sit to stand and pivot transfer to be carried out by a single carer. This reduces the time it takes to carry out the task as well as halving the cost of having two carers present during the transfer. Handle height 790 to 1300mm (31 to 50"). Weight 17kg.

091120781

Maximum user weight



SPARES
Spares are available for this product, for details and full spares listing please contact Customer Services.

1



1



1 Days Orbi-Turn

The Orbi-Turn has been designed for carers to facilitate standing and seat-to-seat transfers, whilst minimising the risk of injury to themselves and the patient. It allows a client with moderate upper and lower limb strength to participate in functional transfers whilst ensuring good working postures for carers and therapists.

- Height adjustable and angled knee pad facilitates normal biomechanical hip and knee movement
- Height adjustable hand bar to suit a range of clients
- Thigh pad gives the client a secure and supportive contact point for increased stability and comfort during transfers
- Turntable effortlessly rotates through 360° for repositioning in even the most confined spaces
- Non-slip footprints provide secure foot placement and show clients where to position feet
- Two rear castors provide easy transport from client to client
- Front mounted sprung castor assists the carer in positioning the Orbi-Turn once the client's feet are in position
- Folds flat for transport and storage
- Height adjustable from 815 to 1315mm to suit a variety of clients
- Steel frame for durability and strength, weighs 15kg

	Lowest Position	Highest Position
Bottom Handle	815mm	1175mm
Top Handle	920mm	1315mm
Knee Pad (Mid Point)	235mm	510mm

091166883

Maximum user weight



SPARES

Spares are available for this product, for details and full spares listing please contact Customer Services.



1 Secure Turn

This new easy-to-use patient turner is ideal for seat to seat transfers when a client's mobility is restricted. It features a curved support frame, which helps to reduce the risk of falling and ensures a safe, comfortable transfer. A hand operated locking device enables the carer to rotate the Secure Turn to any of five positions whilst controlling the turning speed. The non-slip turn disk glides easily in both directions and counter balance points provide stability throughout the transfer. The Secure Turn is fitted with housekeeping wheels so it is easily moved, making it ideal where transfers are required in a number of areas. Width between arms 430mm (17"). Height from handles to platform 760mm (30"). Height of top bar to platform 990mm (39"). Width 660mm (26"). Depth 685mm (27"). Disk diameter 450mm (18"). Weight 17kg.

091169127

Maximum user weight



2 Etac® Molift Raiser Pro

The ergonomic platform for a safe and active sit-to-stand transfer.

Natural movement patterns

- Molift Raiser Pro gives the user the opportunity to use their own muscle power during a transfer. The natural movement pattern is stimulated and the ergonomic design makes the user feel safe during the entire transfer. Molift Raiser Pro is intuitive, which makes it easy to understand and use.

Transfers in confined spaces

- Molift Raiser Pro is stable and can turn within its own radius which makes it very easy to manoeuvre. It is ideal in confined spaces, for instance in toilets. Molift Raiser Pro ensures safe, active standing support during short transfers, such as between beds and wheelchairs/shower chairs. It assists carers when carrying out transfers with minimum risk of injury.

Works in all healthcare environments

- Molift Raiser Pro helps users to train and rehabilitate in all healthcare environments, during every transfer so mobilisation can begin early, right in the intensive care unit. The user can initially improve their balance when sitting and standing up. Molift Raiser Pro always provides safe support.
- Lighter weight - now only 10 kg
- Aluminium base for easier recycling
- New integrated locking brakes with auditive feedback
- Additional fixing point for Raiserbelt & Strap for users of different heights
- New knee support one-handed adjustment with auditive feedback when in position
- Angled knee pad, follows movement of leg when raising from sit to stand
- Hygienic, easy to clean leg support padding (closed cell foam)
- Reduced total width, now 540 mm for narrow doorways
- Lower step on height of only 35 mm
- Height: 1200 mm (47 1/4")
- Depth: 630 mm (25")
- Width: 540 mm (21 1/4")
- Step-on height: 35 mm (1.38")
- Total Weight: 10.7 kg
- Weight: Baseplate 5.4 kg
- Material: Powder coated aluminium, plastic, steel, chrome 3.

091169135

Maximum user weight





1 Medidisc

A small, lightweight stand-on turntable that has a flexible rim allowing the carer to control turning action, giving an effortless transfer. A carrying handle makes this an easily transportable transfer aid. Non-slip studs on underside. Diameter 380mm (15"). Weight 1kg.

091121714

Maximum user weight



2 Transfer Turntable

This strong yet lightweight plastic turntable provides easy rotation through 360° for functional standing and seated transfers.

Ideal for bed to chair or car transfers. Diameter 405mm (15½").

Weight 1kg.

091120831

Maximum user weight



3 Rota Cushion

The Rota Cushion is a specific aid to swivelling patients in a seated position. It consists of two low friction flexible discs that rotate against each other on a central disc. It features an anti-slip base and foam cushioned upper in soft, breathable, waterproof fabric to improve comfort and reduce sweating.

Diameter 380mm (15").

091120666



1 Swivel Seat

This high quality cushioned seat glides you in and out of cars smoothly without struggling. Contours to the seat and will swivel in both directions. May help those who suffer from arthritis, leg or back pain. A strap secures the position of the swivel seat. Diameter 400mm (15½"). Height 70mm (2¾"). Weight 1.6kg. **091121698**

Maximum user weight



2 Revolving Seat

This swivel seat provides invaluable help for getting in and out of the car. The fleece covered foam padded seat is mounted on a wooden ball-bearing turntable that allows the body to rotate through 90°. The cushion has a detachable strap that attaches safely onto the seatbelt. The flame-retardant fleece cover is fully washable. Diameter 380mm (15"). Weight 1.1kg. **091120732**

Maximum user weight



3 Turning Cushion

This foam padded grey cushion is ideal for assisting with seated transfers. It rotates through 360° for smooth and easy moving in any direction. Recommended for chair and car transfers. Comes complete with a removable cover that is machine washable. Weighs 1.1kg and measures 390mm (15¼") in diameter. **091120815**

Maximum user weight



4 Economy Rotary Cushion

A soft padded, flexible turning disc that is ideal for use on car seats, chairs and beds. It allows a smooth swivel transfer without painful jarring or twisting. The underside is covered in an anti-slip material for added safety and the top is cushioned for comfort. Diameter 430mm (17"). **091169267**





1 SuperPole™ with SuperBar™

The SuperPole with SuperBar provides horizontal as well as vertical support, for extra stability. It moves smoothly when lifted and locks securely when let down. The SuperBar locks every 45° and is height adjustable. Supplied as a complete unit with SuperPole: Length 2.44m (96"). Diameter 63mm (2½"). Weight 11kg. Accessory: Length 533mm (21"). Weight 3.2kg.

091121037 Complete Kit

Maximum user weight

136
kg **21 1/4**
st

2 SuperPole™ Basic

The SuperPole Basic is central to a support system providing assistance for standing, transferring or moving in bed. The pole is easily installed using a simple jackscrew expansion and fits ceiling heights from 2.36m (93") to 2.54m (100"). A quality handgrip provides a secure and comfortable grip and rubber pads prevent marking of the ceiling or floor. The SuperPole basic serves as a mounting support for the SuperBar (above) and Safety Island (opposite). Diameter 63mm (2½"). Weight 8.2kg.

081074566

Maximum user weight

136
kg **21 1/4**
st

3 Advantage Rail

This rail has the added feature of a horizontal bar for further support. The locking mechanism on the horizontal bar provides excellent stability while pulling laterally. It can be unlocked by lifting it gently and locked again by lowering it into one of eight positions. Available as a permanently fixed pole or as a portable model that can be stored away or used in another location using the additional floor plate. Height 762 to 965mm (30 to 38"), Handle length 550mm (21"). Weight 10.4kg.

081107440

Maximum user weight

136
kg **21 1/4**
st

1 Homecraft Rope Ladder Bed Hoist

A simple idea to assist with sitting up in bed. The bottom ends of the cords tie around the bed legs, and the plastic ladder rungs enable the user to pull themselves up with a hand over hand movement. Suitable for single or double beds. Rung length 200mm (8"). Diameter 25mm (1"). Cord length 3 metres (118"). Weight 285g.

091089895

Maximum user weight



2 Days Patient Helper

This free standing Patient Helper is placed under the bed base to give a stable overhead support pole for those who find it difficult to rise out of bed. Designed for raising a person to a sitting position, not lifting the entire weight of the person off the bed. The moulded plastic handgrip can be adjusted in height 950 to 1310mm (37 1/2 to 51 1/2"). Disassembles into three pieces for ease of storage and transportation. Minimum required clearance 45mm (1 3/4"). Overall height 1910mm (75"). Overall width 800mm (31 1/2"). Overall depth 610mm (24").

091427400

Maximum load weight



3 Bed Blocks

A pair of lightweight patient hand blocks made from plastic with comfortable angled hand-grips. The high-density foam bases prevent them from slipping on smooth bed covers or causing injury to frail skin on the sides of legs whilst moving. Can be cleaned with a soft brush and detergent in hot water to 75°. Length 215mm (8 1/2"). Width 110mm (4 1/4"). Height 80mm (3"). Weight 610g (Pair).

091121508

4 Swinging Support

Wall mounted support rail assists in horizontal transfers. The support can be set at any of six angles, using a pull chain lock. Strap length adjusts and can also be used in any of four positions along the arm of the support. Length 838mm (33"). Weight 4.7kg.

091024348

Maximum load weight



1



2



3



4



1 Homecraft Mattress Elevator

The Homecraft Mattress Elevator allows those confined to bed for long periods to sit or lie comfortably. It is easily width adjustable for single through to queen size beds. Securely fixed into place using colour coded straps. When in use, it allows the user to rise and recline smoothly, quietly and safely at the touch of a button, thus reducing the need for manual handling. The Homecraft Mattress Elevator allows a user weight of up to 254kg (40st) to be raised and maintained at any angle between 9° and 72°. Available with side guards that, when fitted, will help prevent pets or children from becoming trapped. Overall length 725mm (28½"). Overall width 740 to 1410mm (29 to 55½"). Range of movement 9° to 72°. Weight of elevator unit 11.2kg.

091181668 With Side Guards, UK
091181676 Without Side Guards, UK

Maximum user weight **254 kg** **40 st**



2 Homecraft Bed Grab Rail

An ultra sturdy bed grab rail, offering stability and confidence. The height adjustable frame can be positioned to suit virtually any divan style bed. Two legs rest on the floor and the frame is secured into place via three straps, ensuring minimal movement when in use. The Homecraft Bed Grab Rail is simple to assemble and fit. Height of handle 850 to 950mm (33½ to 37½") from the floor. Width of handle at the top 265mm (10½"). Size of under mattress frame 620x325mm (24½x12¾"). Weight 3.4kg.

091089481

Maximum user weight **127 kg** **20 st**



3 Smart-Rail

Unlike fixed bed rails, the Smart-Rail can be unlocked and pivoted outwards to provide better standing support, meaning less reaching or twisting. The solid and secure steel frame is supplied with straps that are attached to the opposite side of the bed frame to minimise movement during use. Folds flat for transportation and can be easily installed on other beds when travelling. Standard rail height 762mm (30").

081074624 Smart Rail Weight 7.7kg

Maximum user weight **127 kg** **20 st**

4 Parnell Premier Bed Rail

The Parnell bed rail is the ideal solution to help get in or out of bed. Its attractive, ergonomic design complements virtually any bedroom. Functional, strong, hygienic and safe, it provides the user with confidence and independence. Stored flat when not in use, it can be easily assembled and quickly fitted to the bed without the use of tools. Designed for use with electrically adjustable beds, the bed rail can also be used on domestic beds including divans, pine and steel or wooden frame beds.

The bed rail is also compliant with the BS EN 1970:2000 standard. Fits to single or double beds on either side, with a clamp that firmly sits between the floor and the bed base. Comes complete with a full length and width steel floor frame for added stability. The front base is flat so there is no tripping hazard or obstacle to wheelchairs or walking frames. As standard, the rail height is 810mm (31") high. With the optional height extenders, the rail height reaches 925mm (35"). Floor frame length 810mm (31"), width 410mm (16"). Rail height 810mm (31"), width 280mm (11"). Weight 9.5kg.

091172550

Maximum user weight **115 kg** **18 st**

SPARES
 Spares are available for this product, for details and full spares listing please contact Customer Services.

PRICES
 For up-to-date prices please see our Website, alternatively contact your Sales Representative or Customer Services.

1 Uplift Premium Power Seat

The Uplift Premium Power Seat is the portable electric lifting seat that provides 100% lift for those up to 136 kg. For use on most armchairs or sofas, it features uses patented LevelLift™ Technology to gently raise and lower the user.

- Multi-density cushion with memory foam core helps prevent pressure sores
- Zippered, black waterproof cover with comfort mesh for improved airflow, machine washable
- Flexible handle control with large, easy-to-use toggle switch
- Two ergonomic carry handles

Seat depth 497mm (19 ¼"), Seat width 450mm (17"), Height of seat 75mm (3"), Maximum height of seat 355mm (14"), Weight 6kg Cable length 1.2m

091088772 UK plug

Maximum user weight



2 Uplift Seat Assist

Portable, lightweight, lifting cushion that transforms any chair into an automatic lifting chair. Self-powered and easy to adjust, it provides a stable, gentle lift that supports up to 70% of a user's weight. It has a built in carry handle. The memory foam cushion conforms to the shape of the user's body providing a superior level of comfort. Supplied with a black, waterproof cover that can be wiped clean with a cloth or removed and machine washed. Available in two user weight ranges. Seat width 430mm (17"). Seat depth 480mm (19"). Base width 380mm (15"). Base depth 410mm (16"). Seat height closed 50mm (2"). Seat height open 260mm (10"). Weight 3.6kg.

User Weight Range

091088699 35-105kg (5 ½-16 ½st)

091088707 90-160kg (14-25st)

3 UpEasy Lifting Cushion

Similar to the Uplift Seat Assist but with a polyurethane foam cushion and a navy, washable cover. Dimensions are the same as the Uplift Seat Assist.

User Weight Range

091088673 43-100kg (6 ¾-15 ¾st)

091088681 90-154kg (14-24 ¼st)

4 Portable Raised Seat

This seat is helpful when transferring to and from a sitting position. The lightweight polyurethane/polystyrene foam cushion has a carry strap, making it easy to lift and carry. Comes with removable and washable cotton cover. Measures 375×375mm (14¾×14¾"). Latex free.

Seat Height


081174903 100mm (4")

Maximum user weight





days

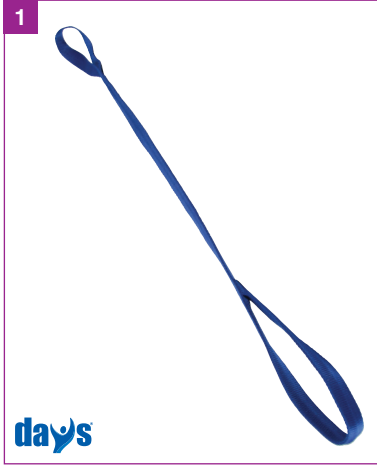

days



days



days



1 Days Leg Lifter

This webbing strap with stiffened stem enables the user to move a stiff or immobile leg when using a bed, wheelchair, footstool etc. The upper loop fits around the hand and wrist, and the lower stiffened loop slides over the foot to give control in lifting the leg. Length 864mm (34"). Weight 78g.

091118793



2 Days Leg Lifter with Stirrup Loop

Similar in style and function to the leg lifter above, but with a stiffened triangular stirrup, making it easier to fit over shoes and giving more support to the foot. Length 915mm (36"). Weight 120g.

091118801

3 Days Dual Handle Leg Lifter

Similar in construction to the Leg Lifter (091118793) above, but with two handles for users that cannot lift their leg with one hand. Length 865mm (34"). Weight 110g.

091118819

4 Dycecm™ Non-Slip Floor Mat

Dycecm's nonslip floor mats are perfect for providing a sure footing, boosting mobility and helping with any balance issues. These mats are manufactured using Dycecm's patented non-slip material, meaning that you won't get the unbeatable grip and stability anywhere else. Dycecm's floor mats can help in lowering the risk of falls, and have been used by therapists to improve standing, sitting and positioning. Gain a sure footing with Dycecm's non-slip floor mats. Length 600mm (23½"). Width 450mm (18"). Weight 943g.

091102912 Blue

091102920 Grey

5 Ambulance Chair

This ambulance chair features front and back carrying handles and can be folded flat when not in use. Manufactured from plastic coated steel it has hygienic wipe clean polyester upholstery, a nylon webbing safety strap and non-marking wheels. Handle Height 74cm (29"), Seat Size 35x29cm (14" x 11½").

Weight 9kg.

091439983

Maximum user weight

102 kg **16½** st

6 HandyBar®

Simple, lightweight stainless steel bar with soft, non-slip hand grip, to provide strong, secure assistance in getting in and out of a car. Acts as a grab bar without the need for costly or unsightly adaptations. Only 300g and 230mm (9¼") long, it is easy to store in the car or carry along. Fits the majority of vehicles with a u-shaped striker.

091032705

RETAIL PACKAGED
These products are supplied in retail packaging.